

## MINIWALK®



*MINIWALK® size no. 1.*



*MINIWALK® size no. 2 with accessories.*

MINIWALK® from Meyland-Smith in Denmark is the newest model of walking supports for small children from approx. 1½ to 7 years of age. The walking support is based on the famous lines and qualities of MEYWALK®. The height of the trunk support and seat are adjusted quickly and stepless with spanner grips. MINIWALK® is available in two sizes, where size no. 1 is for children from 1½ to 4 years of age and size no. 2 is for children from 4 to 7 years of age. The seat and trunk support has spring suspension to encourage a better and more natural gait. The body posture of the child can be set in any position from

vertical to forwards tilt by 15°, making it possible to find the body inclination that gives the child the best control. The trunk support is obtainable in two different sizes, and the springs are obtainable in three different strengths. MINIWALK® is made with materials of high quality, that guarantee long life, maximum reliability and a minimum of maintenance. MINIWALK® is therefore suitable for indoor as well as outdoor use.

MINIWALK® comes complete with seat, trunk support, rear stop, handle bar, fender wheels and parking brakes.

By replacing the two uprights (the tubes with the rulers) MINIWALK® can be changed from size no. 1 into size no. 2. In this way MINIWALK® grows with the child.

The seat and trunk support on MINIWALK® has spring suspension, and the springs are obtainable in three different strengths, which should be **selected according to the weight of the child:**

- Min. 10 kg (22 lbs.) = Green springs
- Min. 15 kg (33 lbs.) = Silver springs
- Min. 20 kg (44 lbs.) = Black springs

**Please inform spring strength when ordering!**

The trunk support on MINIWALK® is obtainable in two different sizes, which should be **selected according to the measure of the child's thorax:**

- Max. 70 cm (27") in circumference
- Max. 90 cm (35") in circumference

**Please inform trunk support size when ordering!**

## Construction

- Frame manufactured from steel tubes, powder coated.
- Components and accessories chromium plated / electro-galvanised.
- Solid (purbane) wheels.
- Rims and front wheel forks in black hard plastic with ball bearings.
- Paddings of black water-resistant foam rubber coating.

(We reserve the right to make changes in design and construction.)

## Technical data

MINIWALK®	Size no. 1		Size no. 2	
Length:	74 cm	(29")	74 cm	(29")
Width:	65 cm	(25,5")	65 cm	(25,5")
Seat/groin height:	27 - 48 cm	(10,5" - 19")	27 - 63 cm	(10,5" - 25")
Trunk support height:	52 - 75 cm	(20,5" - 29")	67 - 90 cm	(14,5" - 35,5")
Wheel size, front:	150 x 30 mm	(6" x 1¼")	150 x 30 mm	(6" x 1¼")
Wheel size, rear:	250 x 50 mm	(10" x 2")	250 x 50 mm	(10" x 2")
Weight:	15 kg	(33 lbs.)	15 kg	(33 lbs.)
Max. load:	40 kg	(90 lbs.)	40 kg	(90 lbs.)

## Accessories


- Anti-tip supports
- Drag brakes
- Non-reverse brakes (can be deactivated)
- Swivel locks for front casters (can be deactivated)
- Leg guides
- Hip pads
- Crossbar - no seat option
- Leg separation plate



**MEYLAND-SMITH A/S**

INDUSTRIVEJ 27  
DK-9830 TAARS

Tel: +45 98 96 19 85 • Fax: +45 98 96 19 86  
www.meyland-smith.com • info@meyland-smith.com

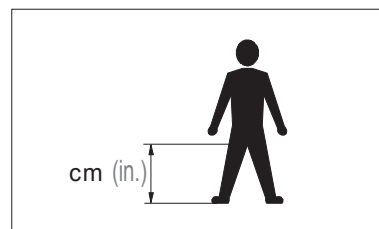

*This is how the seat/groin height is measured.*

## Frame colours



*Turquoise (RAL 5021)*



*Red (RAL 3002)*

## Padding colour



*Black*



*Please get in !*

